



INNOVATING DIABETES CARE




Why Equil?

Tubeless
Leave your diabetes behind

Reliable
Stay secure with our multiple alarms

Compact
23g, as light as 5 nickels


Intelligent
Freedom of insulin delivery anytime, anywhere

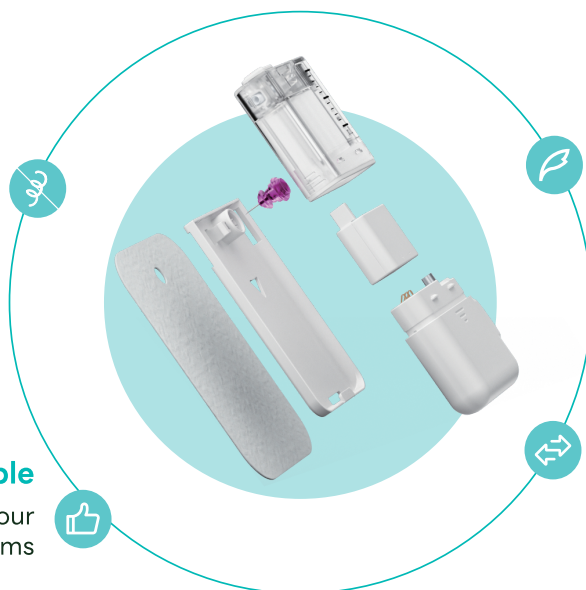

Pump function


Mimicking the insulin release of a healthy pancreas.

Mode 1: Bolus

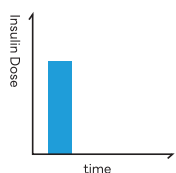


Meet daily meal needs

Mode 2: Extended Bolus (Dual-wave bolus)

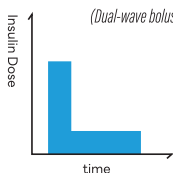


When you eat meals with both rapidly and slowly absorbed carbohydrates

Mode 3: Square Bolus (Square-wave bolus)

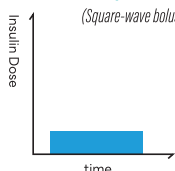

When You eat meals that last for a long time

Mode 4: Temporary Basal

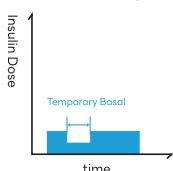

Meet the need for activities



► Bolus Calculator

The Bolus Calculator assists users in determining the appropriate bolus size considering factors such as blood glucose levels and other relevant parameters.

Extended Boluses

Utilize an Extended Bolus with Amount Now for meals containing a mix of rapidly and slowly absorbed carbohydrates, including those with high-fat content like certain fast foods. This bolus is delivered in two stages – an immediate amount followed by the remainder over an extended period.

Temporary Basal

Typically, basal rates remain stable, but under specific circumstances, adjusting to a temporary basal rate may be necessary to prevent fluctuations in blood glucose levels.

Specification



Basal Rate
0-35U/h



Basal Programs
7 Basal Programs, each with 48 time segments



Maximum basal rate
0.1-35U/h, Default: 1.5U/h



Quick bolus
On/Off, Default is on



Base basal rate
0-35U/h, Default: 0.5U/h



Temporary basal rate
U/h or % of basal rate, last 3 basal rates remembered, Default is Off



Bolus size
0.025U-35U, 3 Presets



Quick bolus increment
0.1-2U



Bolus increment
0.025/0.05/0.1/0.5/1U, Default: 0.1U




Max bolus size
1-35U



Extended bolus
Programmed in U or % of total bolus, Default is off, Extend time: 0.5-8 hours in increments of 0.5 h

Operation Steps

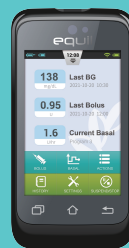

PDA



Pump



Cannula Inserter



Cannula Pump base




Reservoir

Component

Consumables


1.

Fill the reservoir

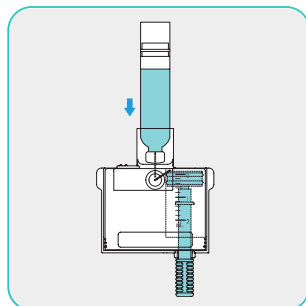


2.

Apply the infusion set

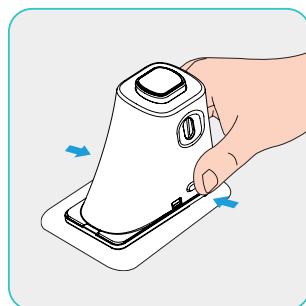


3.

Attach the pump

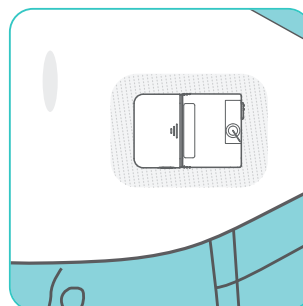

4.

Control with PDA

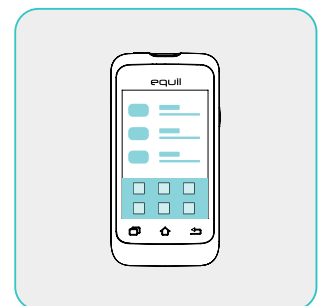

Manufactured by:

MicroTech Medical (Hangzhou) Co., Ltd.

No.108 Liuze St., Yuhang District, Hangzhou, 311121 Zhejiang, P.R.China

✉ globalmarketing@microtechmd.com

🌐 www.microtechmd.com/en/



Local Authorized Representative:

Nano Medic Care Sdn. Bhd.

📍 1-12, 1st Floor, Jalan Indah 2, Taman Indah, 84000 Muar, Johor, Malaysia

☎ Phone: +606-9513111 | 📠 Fax: +606-9546111 | 📅 Monday-Friday: 9am-5pm

📞 Customer Service Hotline: +6018-9704732



SCAN HERE FOR MORE INFO!